


## Plotting example snippets of EEG data

This script loads a single recording file using the MNE package (v1.6.1) in python (v3.12) and plots a 5-10 second snippet of data from 4 selected channels using Matplotlib (v3.8.3). The output of the script should appear as the figure below saved as a .png into the output folder specified in the script.

Output file:

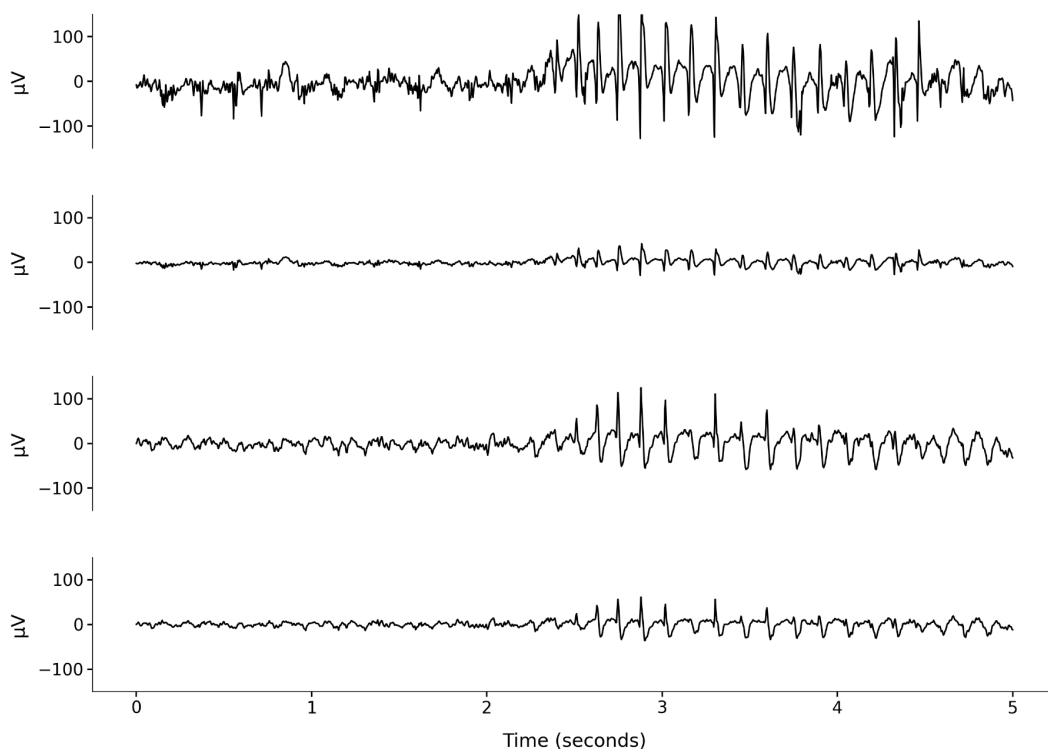

## Instructions:

1. Navigate to the main function and enter the recording file you want to use in the recording\_folder parameter

```
recording_folder =  
'/Volumes/Sarah/SYNGAPE8/DATA/SYNGAPE8/12W/SYNGAPE8_2780/TAINI_1048_2780_EM4-  
2024_04_05-0000.dat'
```

2. Edit the output path to the desired folder you want to save the output figure into

```
output_figure_path =  
'/Volumes/Sarah/SYNGAPE8/OUTPUT/SYNGAPE8/12W/SYNGAPE8_2780/'
```

3. Enter the EEG channels you want to plot in `eeg_channels` parameter

```
eeg_channels = [5, 8, 9, 10] # change these to the channels you want to plot
```

4. Navigate to the `crop_by_start_and_end` function and enter the start and end time in seconds of the snippet of EEG data you want to plot. `tmin` is the start or minimum time and `tmax` is the end of max time.

```
cropped_raw = custom_raw.crop(tmin=137943, tmax=137947) # crop up to 10  
seconds for example plotting
```

Note : if you are not sure about the section of EEG data you want to plot, plot the whole recording in an interactive plot and search for a snippet you want to plot and note the start and end times in seconds

5. If you want to edit the name of the .png saved, navigate to the `plot_example` function and enter in the name you want the file to be saved as

```
plt.savefig(output_figure_path + '/' + 'ExampleTrace_0100.png', dpi=200)
```